

# Classification Of Arrhythmia By Using Deep Learning With 2-D ECG Spectral Image Representation

## MODEL BUILDING

### INITIALIZING THE MODEL

Team ID	PNT2022TMID31222
Project Name	Classification Of Arrhythmia By Using Deep Learning With 2-D ECG Spectral Image Representation

### INITIALILIZING THE MODEL:

keras has 2 ways to define a neural network:

- Sequential
- Function API

The Sequential class is used to define linear initializations of network layers which then, collectively, constitute a model. In our example below, we will use the Sequential constructor to create a model, which will then have layers added to it using the add () method.

Now, will initialize our model.

11/7/22, 12:35 AM

Untitled8.ipynb - Colaboratory

#### ▼ Importing Keras libraries

```
import keras
```

#### ▼ Importing ImageDataGenerator from Keras

```
from keras.preprocessing.image import ImageDataGenerator
```

#### ▼ Importing Keras libraries

```
[1] import keras
```

#### ▼ Importing ImageDataGenerator from Keras

```
[13] from matplotlib import pyplot as plt
      from keras.preprocessing.image import ImageDataGenerator
```


#### ▼ Defining the Parameters

```
train_datagen=ImageDataGenerator(rescale=1./255, shear_range=0.2, rotation_range=180, zoom_range=0.2, horizontal_flip=True)
test_datagen=ImageDataGenerator(rescale=1./255)
```

```
<keras.preprocessing.image.ImageDataGenerator at 0x7fb7448ac110>
```

## APPLYING ImageDataGenerator to train dataset:

`flow_from_directory ( )` method for Train folder.

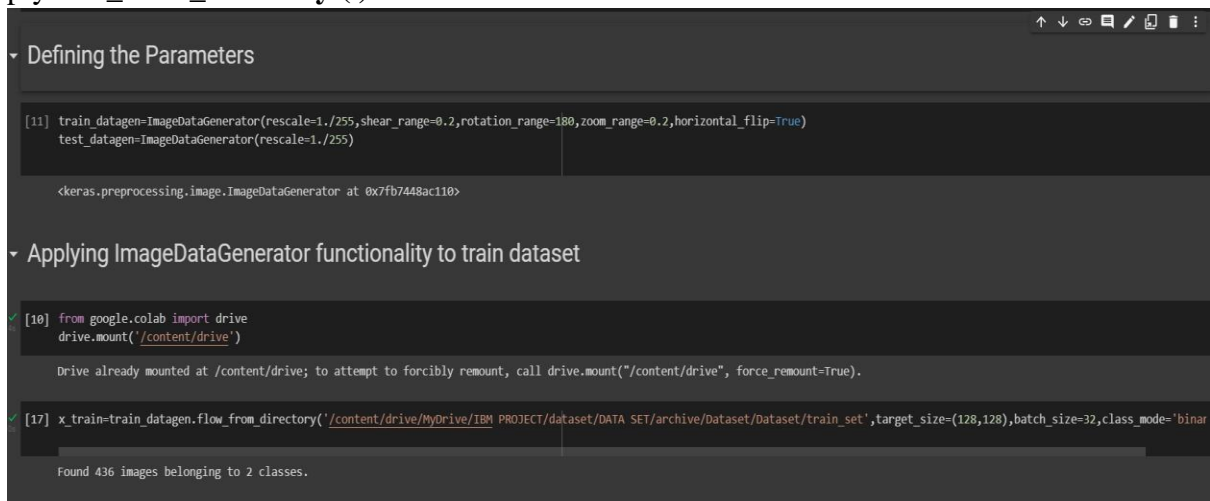

```

# Defining the Parameters

[11] train_datagen=ImageDataGenerator(rescale=1./255, shear_range=0.2, rotation_range=180, zoom_range=0.2, horizontal_flip=True)
     test_datagen=ImageDataGenerator(rescale=1./255)

<keras.preprocessing.image.ImageDataGenerator at 0x7fb7448ac110>

# Applying ImageDataGenerator functionality to train dataset

[10] from google.colab import drive
     drive.mount('/content/drive')


Drive already mounted at /content/drive; to attempt to forcibly remount, call drive.mount("/content/drive", force_remount=True).

[17] x_train=train_datagen.flow_from_directory('/content/drive/MyDrive/IBM PROJECT/dataset/DATA SET/archive/Dataset/Dataset/train_set', target_size=(128,128), batch_size=32, class_mode='binary')

Found 436 images belonging to 2 classes.
```

## APPLYING ImageDataGenerator to test dataset:

Applying the `flow_from_directory ( )` method for test folder.

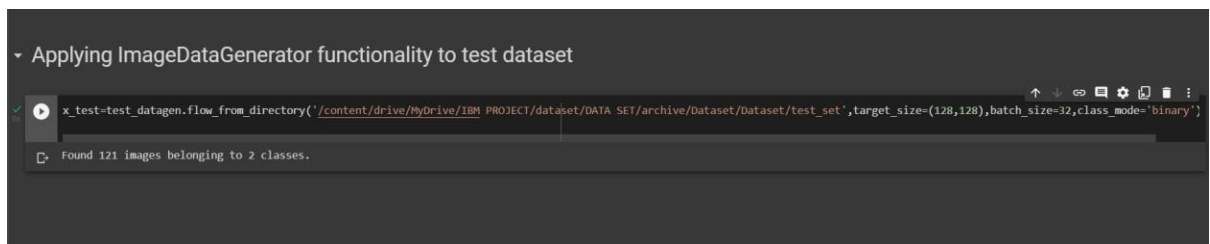

```

# Applying ImageDataGenerator functionality to test dataset

x_test=test_datagen.flow_from_directory('/content/drive/MyDrive/IBM PROJECT/dataset/DATA SET/archive/Dataset/Dataset/test_set', target_size=(128,128), batch_size=32, class_mode='binary')

Found 121 images belonging to 2 classes.
```

## IMPORTING MODEL BUILDING LIBRARIES:

11/8/22, 1:16 AM

Main code - Colaboratory

### ▼ Importing Model Building Libraries

```
#to define the linear Initialisation import sequential
from keras.models import Sequential
#to add layers import Dense
from keras.layers import Dense
#to create Convolutional kernel import convolution2D
from keras.layers import Convolution2D
#import Maxpooling layer
from keras.layers import MaxPooling2D
#import flatten layer
from keras.layers import Flatten
import warnings
warnings.filterwarnings('ignore')
```

## INITIALIZING THE MODEL:

### ▼ Initializing the model

```
model=Sequential()
```